GB

Acupressure wristbands

FR

Bracelet d'Accupression Anti-nausée

ES

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NL

Acupressuurpolsbandje

DE

Akupressur-Armband

IT

Bracciale per digitopressione

PL

Opaska na nadgarstek do akupresury

SE


Akupressur armband

CZ

Akupresurní náramek

JP

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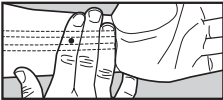Apothecary Products, LLC

11750 12th Ave S., Burnsville, MN 55337, USA

 (800) 328-2742  info@apothecaryproducts.com

ApothecaryProducts.com

Imported by: SANICO S.r.l.

Via G.Ferraris 37/41 - Cusago (MI) - ITALY

 info@sanicare.it - www.sanicare.it

CE

EC

REP

MD

MADE IN CHINA

EMERGO EUROPE

Westervoorstsedijk 60

6827 AT, Arnhem

The Netherlands

04032-A

GB

ENGLISH

Product Description: Acupressure wristbands

Intended Use: Acupressure wristbands are intended to stimulate the P6 Neiguan acupressure point and help in the relief of nausea and vomiting caused by motion sickness, anesthesia or chemotherapy.

Instructions / Directions: The Acupressure wristbands operate by exerting pressure on the P6 Neiguan (acupressure) point on each wrist. To find your P6 Neiguan point, place your middle three fingers on the inside of your wrist with the edge of the first finger just below the first wrist crease. The P6 Neiguan point is located underneath your third finger in between the two wrist (flexor) tendons. Position the button facing inward over the P6 Neiguan point. Repeat the process for the other wrist. One band must be worn on each wrist to be effective. In case of children locate the point with their own fingers. Effectiveness may vary with individual patient characteristics. It may take up to five minutes to become effective, but usually lasts as long as you wear the bands. For better results put bands on before starting travel. Greater relief can be obtained by applying finger pressure directly to the button. The Acupressure wristbands should only be used on the P6 Neiguan point. Use at any other location will not relieve nausea. The Acupressure wristbands may be washed with mild detergent in warm water up to five times. Replace when the bands no longer fit snug around the wrist.

Precautions/Cautions: Discontinue use if product causes discomfort. Do not use if broken or damaged. The Acupressure wristbands provide symptomatic relief, it will not cure or remove the cause of the nausea. Nausea and vomiting can be associated with serious conditions. If your symptoms persist, consult a physician at once. If skin irritation occurs beneath the band, or if pain or swelling occurs, remove the bands immediately and consult a physician if these symptoms persist. Store out of reach of children. The button is a choking hazard if removed. Individuals with known or suspected circulation problems should consult a physician prior to use. If, in relation to the use of Acupressure wristbands, a death or a serious deterioration of health has occurred, this should be reported to the manufacturer and the competent authority of your country.

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