
SANICO

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MATERASSO ANTIDECUBITO AD ARIA
ANTI-DECUBITUS AIR MATTRESS
ANTI-DECUBITUS LUCHTMATRAS
MATELAS PNEUMATIQUE ANTI-ESCARRES
ANTI-DEKUBITUS-LUFTMATRATZE

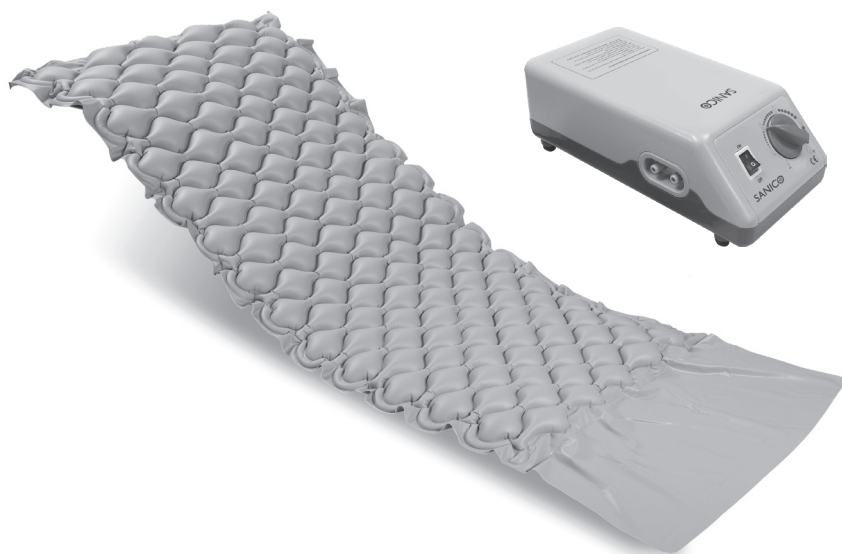IT EN NL FR DE

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USER'S MANUAL

Anti-bedsore alternating pressure bubble mattress

1. Introduction
2. Product description
3. Application
4. Contraindication
5. Use
6. Safety of use **WARNING!**
7. Everyday use and maintenance
8. Technical parameters
9. Symbols
10. Contact details

PLEASE READ WHOLE INSTRUCTION CAREFULLY BEFORE USE.

1. INTRODUCTION

Thank you for choosing our product. We are deeply convinced that it will meet your expectations. Please read this document before first use. This manual is passed in order to provide the necessary information for the safe use of the device. The correct application of the recommendations will contribute significantly to extend the durability and aesthetics of the product. We want to emphasise that constantly we are trying to improve the quality of our products, which may result in minor changes that are not included in the manual.

2. Product description

Items that are included in the package:

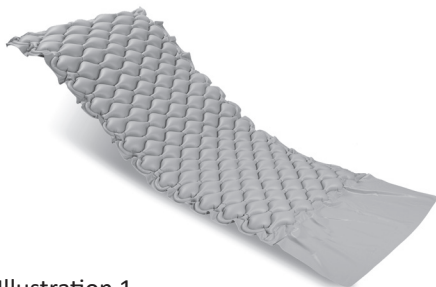
Illustration 1

Illustration 2

Illustration 3

Alternating pressure bubble anti-bedsore mattress equipped with end flaps which allows to place the mattress on the bed. Maximum weight capacity: 130 kg

3. Application

Prevention and support of bed sore treatment

4. Contraindication

In case of adverse skin changes you should immediately consult a physician or medical personnel.

5. Use

- Place the mattress on a standard bed so that the air intake would be located on the side of patient's legs. Put the end flaps under the bed mattress.
- Check if the switch is in the "off" position.
- Connect the pump to electricity

Illustration 4

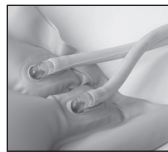
Illustration 5

- Connect the pipes to the pump – insert them into the air outlet located on the pump side (Illustration 4.) and on the mattress (Illustration 5.).
- Place a pump on hard surface or hang it on hinges. The place where the pump will stand must be clean and dry. Dust, hair or other impurity may block the filter.
- Please, check if the pipes are not folded or pressed.
- Switch on the pump.
- Set the pressure regulator on maximum. It will shorten the time of pumping which usually takes 20-30 min.
- After inflating the mattress, you should choose the pressure according to the doctor's recommendation, taking into account that generally the pressure of the mattress should be proportional to the weight of the patient, also considering his comfort.
- Place a patient on the mattress.

In order to receive exact instructions concerning an appropriate use and optimal pressure, please, consult your doctor. The mattress and pump create effect of "constant movement" (changing position of body support points), provide constant muscle massage and improve blood circulation. Thus, the mattress is an effective tool to prevent bedsores. But the position of the patient on the mattress needs to be changed despite of using it.

Switching off

- Turn the switch into the "OFF" position.
- Disconnect the pipes from the air outlet in order to prevent air from escaping.

6. Safety of use WARNING!

- a) Before using the mattress, refer to this manual and keep it for future reference.
- b) Permitted user weight 130 kg.
- c) Please remember that on the mattress should stay only the patient (other people on the mattress are not allowed). It may cause permanent damage to the mattress.
- d) You are not allowed to do the physiotherapy treatment or rehabilitation on inflated mattress.
- e) The mattress should be protected against contact with body care cosmetics. It may cause permanent damage of the mattress. It is recommended to use an insulating sheet, that is, one that does not let liquids through. There are no sheets in the mattress equipment.
- f) You must not use the mattress close to inflammable materials. It can lead to explosion!
- g) Always check if the pump is plugged into appropriate socket (110-220 V/50 Hz).
- h) The plug can not be removed from the socket when the pump is on.
- i) Always make sure the pump is in a dry place, remember that it cannot get wet.
- j) Do not touch the pump when it is wet. Always unplug the pump from electricity when it is wet.
- k) Never leave working pump unattended.
- l) Pay special attention if the pump is used in the presence of children or people with disabilities.
- m) Pump can be used according to its purpose only.
- n) Do not use it with other, unauthorized devices.
- o) Never use the pump with damaged or moist cable / plug. You should submit such a pump to a service point.
- p) The cable should be kept away from the hot surfaces.
- q) Never leave the pump in a place where the cable or tubes can be dangerous for people who are not familiar with the device.
- r) Always pull the plug not the cable while disconnecting.
- s) Never block the pump filter. The filter must be free from any impurities such as hair, scarfskin etc.
- t) Always consult your doctor before using the mattress.
- u) Using the isopropyl alcohol or dissolvent may damage the pump.
- v) Please, make sure that the pump is completely dry before setting on.
- w) The mattress can be folded while storing. Before folding, check if it is completely dry which protects it against mold. The mattress should be stored away from light and in a dry place.
- x) Do not put the mattress close to the fire.
- y) Keep the mattress away from sharp objects.
- z) Don't put the mattress directly on springs or frame with sharp edges.
- aa) Pay attention to skin changes. If necessary, consult your doctor or nurse.
- bb) Tubes cannot be bent.
- cc) Never use the mattress deflated.
- dd) While moving bed, unplug the mattress.
- ee) The mattress can be used on the flat surfaces only. Using the mattress on the orthopaedic adjustable beds it is not recommended.

7. Everyday use and maintenance

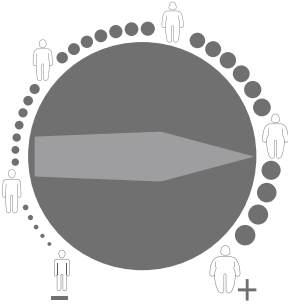
Illustration 6

L'illustration 6 shows the range of adjustment mattress inflating located on the pump. Depending on the doctor's indications, your weight and your preferences, you can choose the right degree of firmness by acting on the knob.

The mattress is inflated by a system of pipes. The cells are alternately inflated which ensures changing position of body support points. Patient's body rests on the same, invariable height regardless of the phase of pump functioning

Phases of pump functioning:

- 1 phase: Inflating all the cells
 - 2 phase: Deflating the first section of cells
 - 3 phase: Reinflating all the cells
 - 4 phase: Deflating the second section of cells
- The phases occur in cycles every 12 minutes.

PROBLEM SOLVING

Do not try to repair the pump, do not open it as it can lead to electric shock!

The knob is in "ON" position but the pump does not work.

- Please, check if the cable is plugged to electric supply and to the pump. If not, you should plug it.
- Please, check if there is a voltage in the electric circuit. If there is not, find the cause of such situation.
- Please, check if the cable is not damaged. If so, please contact the manufacturer.
- The pump works too slowly or the pressure is too low.
- Please, check if the pressure regulator is set on an appropriate pressure. If not, please, regulate it.
- Please, check if the pipes are correctly connected to the air outlet. If not, press them in such a way which prevents the air from escaping.
- Please, check if the pipes are not folded or pressed. If so, please, unbend them.
- Please, check if the cells are not damaged. If so, please, repair it using the reparation kit.
- Please, check if the pump works correctly – if you can hear some strange sounds like hissing etc., please contact the manufacturer

CARE AND MAINTENANCE

Disconnect the pump from electricity. Clean all the surfaces with a moist cloth. Before switching on it has to be dry.

The mattress should be cleaned with water with soap. After cleaning, dry it carefully. The mattress mustn't be washed in high temperature, dried in a washing machine, bleached, chemically washed, ironed or pressed.

The mattress may be disinfected using a product Incidin Plus, in accordance with the instructions for use of the preparation.

You must not use:

- solvents,
- toilet detergents,
- sharp brushes and hard objects,
- detergents that contain chlorine, corrosive detergents,
- Do not clean the exerciser with pressure or steam aggregates!

Frequency	Object of control	Comments
Before first installation	Checking the technical condition and all the elements within the packaging.	Check if any elements have any visible mechanical damage.
Once a month	Cleanliness and general condition.	In case of mechanical damage contact the producer immediately.

Repair

In case of reparation, please contact the producer.

Customer service

If you have any questions or need help, contact the producer.

8. Technical parametres

Pump:

Power supply: 110/220V 50Hz

Air output: 5-7 litri/min.

Pressure range: 50 mmHg - 120 mm Hg

Cycle time: 12 min.

Power consumption: 7 W

Size: 23,5 cm x 13 cm x 9,5 cm

weight: 0,8 kg

Pump volume: 30-35 dB

Mattress:

Size: 200 cm x 90 cm x 70 cm

weight: 1,8-2,1 kg W

eight: 130 Kg

Series number: on the packaging

Durability: 3 years from production date

9. Simboli

Confirmation of compliance with EU standards

Manufacturer

Date of production

Read the user's manual

Warning

Catalog number

Serial number

Protect from humidity

Protect from sunlight

Permitted user weight

Use inside buildings

The product is intended to transport by car

The product is intended to the transport by plane

Medical device

DISPOSAL OF YOUR OLD APPLIANCE

When this crossed-out wheeled bin symbol is attached to a product, it means the product is covered by the European Directive 2002/96/EC.

All electrical and electronic products should be disposed off separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.

The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

CONDITIONS FOR SERVICE UNDER GUARANTEE

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 2 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: electrode cups, elastic bands, power supply cables etc. nor regular checks and maintenance (calibration).
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@sanicare.it.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases:

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;
- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect use, installation or transport of the product;
- components or accessories subject to rapid wear and tear;
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and profit-making purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

e-mail: info@sanicare.it

www.sanicare.it

Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

Send this coupon in case of repaires:

Product: PL023

Serial Number:

if available

Date of purchase:

		/			/				
-------------	-------------	---	-------------	-------------	---	-------------	-------------	-------------	-------------

Dealer's Stamp:

--

Buyer's Full Name:

Street/Square:

N°:

City and State:

Postal Code:

Country:

Phone Number:

E-mail:

@

Problem description:

Signature:

Date:

		/			/				
-------------	-------------	---	-------------	-------------	---	-------------	-------------	-------------	-------------

WARRANTY IS VALID ONLY IF ACCOMPANIED BY INVOICE/TICKET.

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Foshan Hongfeng Co., Ltd
MedNet EC-REP GmbH
SANICO Srl - Via G.Ferraris, 37/41 - 20047 Cusago (MI) ITALY