SANICO

AEROSOL A PISTONE
COMPRESSOR NEBULIZER
NÉBULISEUR DE COMPRESSEUR
KOMPRESSOR-NEBULIZER
COMPRESSOR VERNEVELAAR
NEBULIZADOR DE PISTÓN

PL044

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Indice

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Table of Contents

1. IEC Symbols
2. Important Safeguards
3. Introduction
4. Specifications
5. Description
6. Accessories
7. Operation Instruction
8. Cleaning
9. Maintenance
10. Technical data
11. EMC
12. Troubleshooting
13. Disposal
14. Manufacturer
15. Warranty

IEC Symbols

Caution Attention or Consult accompanying documents

Type BF Equipment

Alternating Current

Double Insulation

Indoor use only

Important Safeguards

Caution:

1. Follow the instruction of your physician to operate this unit.
2. This product is a nebulizer for the inhalation of medical aerosols and is suitable for children used solutions. Use only the type and amount of medication prescribed by patient's doctor.
3. This product is intended for aerosol therapy only. And other use is not recommended.
4. Do not untwist the jet nebulizer during operation.
5. Do not pour more than 6ml of solution into the jet nebulizer.

Danger: To Reduce the risk of electrocution.

1. Always unplug the unit immediately after using.
2. Do not use while bathing.
3. Do not place or store the unit where it can fall or be pulled into a tub or sink.
4. Do not place or drop into water or other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately.

Warning :To Reduce the risk of burns, electrocution, fire or injury

1. Electrical Shock Hazard—do not remove the cabinet or open the cover.
2. Disconnect the power cord from the electrical outlet before cleaning or servicing.
3. Do not place this equipment near hot, sparking or burning objects.
4. Do not use oil or grease on or near this device.
5. Turn the unit off when it is not in use.
6. Keep the cord away from HEATED or HOT surfaces.
7. NEVER drop or insert any object into any opening.
8. NEVER block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked.
9. Avoid operating in wet or damp locations.
10. Unplug the unit before filling the nebulizer.
11. When using this product near TV, microwave oven ,pulse telephone X ray or other strong electric field, it will be disturbed, suggested to be far away from these device measuring.

Warning: To reduce the risk of infection:

1. This Nebulizer kit is intended for both adult and pediatric patient use.
2. Cleaning of the nebulizer is recommended after each aerosol treatment. disinfecting is recommended once a day. Please follow cleaning and disinfecting instructions in this manual.

This unit is not suitable for use in anaesthetic breathing systems

This unit is not suitable in suspension or high viscosity form. In such cases information should be sought from the drug supplier

Introduction

This compressor nebulizer is designed to deliver the prescribed medication solution to treat patient respiratory disorders, such as asthma, allergies and bronchitis. The nebulizer converts the medication solution into an aerosolized mist which is inhaled by the patient through the mouthpiece or masks. Read this manual thoroughly before using the nebulizer and save for future reference

Specification

Electrical Requirements	230 VAC / 50Hz
Power Consumption	≤160VA
Rated Current	0.70A
Dimensions	Length 145mm × Width 142mm × Height 90mm
Weight	1.16kg
Sound Pressure Level	≤58 dB(A) (at 1meter distance)
Max.Comp.Pressure	≥35 Psi (241 kPa)
Neb.Operating.Pressure	9~16Psi
Liter Flow Range	5-8 L/min
Nebulization Rate:	≥0.18 ml / min (by weight loss)
Particle size	MMAD approximately $1.99 \pm 0.69 \mu\text{m}$
Maximum nebulizer solution capacity	6ml
Maximum residual medicine volume	≤1.0 ml
Operating Temperature range	+5° to +40°
Operating Humidity range	30% R.H. to 85% R.H.
Storage Temperature range	-20° to +60°
Storage Humidity range	10% R.H. to 93% R.H.
Operating atmospheric pressure range	700~1060hPa

Protection against electric shock

-classified as Type II

-Type BF applied parts: Mouthpiece, masks

Degree of safety in the presence of flammable anesthetics or oxygen


—No AP/APG (not suitable for use in the presence of flammable anesthetics or oxygen)

Mode of operation—Continuous

IP21—Protection against harmful ingress of water is Ordinary

Description

- | | | |
|-------------------------|-------------------|-------------------------|
| 1. Nebulizer cup | 4. Pediatric mask | 7. Filter & cover |
| 2. Mouthpiece/Nosepiece | 5. Power switch | 8. Air outlet connector |
| 3. Adult mask | 6. Air tube | |

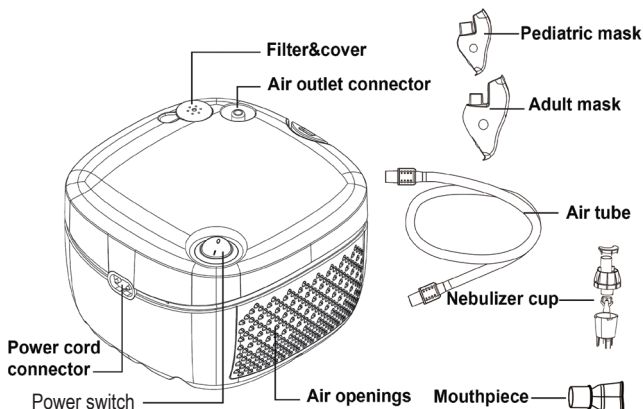Accessories

Nebulizer machine	
Air tube	
Adult mask	
Pediatric mask	
Nebulizer cup	
Mouthpiece/ Nosepiece	
5pcs air filters	

Features and cautions:

The nebulizer cup is adjustable, when the adjusting port opened, the nebulizer rate is fast, when the adjusting port closed, the nebulizer rate is slow. For inhaled medication to be most effective, we suggest the adjusting port closed for pediatric treatment, the adjusting port is opened for adult treatment.

Operation Instructions

1. Place compressor on a stable, sturdy and flat surface so that the unit can be easily reached when you are seated.
2. Pull out the power cord from the base of the unit, and make sure that the unit is in the "off"(O) position by pressing on the left side of switch.
3. Plug power cord into a wall outlet.
4. Connect one end of the nebulizer tubing into the air outlet connector.
5. Add the prescribed solution through the opening on cap using an eye dropper or pre-measured dose container.
6. And the marking "6cc" marked on the nebulizer cup is the maximum level
7. And the medication can not be filled beyond the maximum level
8. Assemble mouthpiece and T-mouthpiece, and insert into the top of the nebulizer cap. If an aerosol mask is used, connect the bottom of mask to the top of nebulizer cap.
9. Connect tubing to nebulizer air inlet connector.
10. Turn the power on to start the unit.
11. Implement the treatment by placing the mouthpiece between the teeth. Inhale and exhale through the mouthpiece.
12. If an aerosol mask is used, place the mask over the mouth and nose.

Cleaning

1. Turn the power off and unplug from the wall outlet .
2. Remove the tubing from the air inlet.
3. TO CLEAN: Disassemble mouthpiece, mouthpiece, nebulizer cup, chamber and baffle, and wash this items in hot water with dish-washing detergent. Rinse these items thoroughly to remove the detergent and let air dry.
4. TO DISINFECT: Mix one part of white vinegar with 3 parts hot water in a clean container. Submerge mouthpiece, mouthpiece, nebulizer cap, chamber and baffle for half an hour in the solution. Remove from solution and let air dry.
5. There is no need to clean the tubing. If necessary wipe the surface regularly.

CAUTION:The jet nebulizer must be changed if it gets clogged.

CAUTION:The nebulizer and accessories must not be boiled.

Maintenance

1. The filter of air compressor should be replaced when it turns gray. Additional filters can be procured from your dealer.
2. Only authorized personnel shall make repairs to this product.

CAUTION: The nebulizer may be damaged if used with a dirty filter or if the filter is replaced by other material such as cotton. Do NOT operate without a filter.

Technical data

Particle Size: MMAD* approximately $1.99 \pm 0.69 \mu\text{m}$


Appropriate Medication Quantities: 2ml minimum - 6ml maximum

Aerosol Output: $4.48 \pm 0.32 \text{ ml}$

Aerosol Output Rate: 0.06ml/min (2ml, 1%NaF)

Compressor nebulizer the particle size distribution curve as well:

Result of cascade impactor **measurements for particle size with the Compressor nebulizer and nebulizer kits

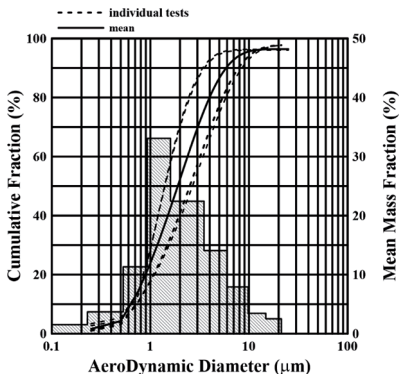*MMAD = Mass Median Aerodynamic Diameter

Important information regarding Electro Magnetic Compatibility (EMC)

1* WARNING: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.”

2* WARNING: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.”

3* WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the PL044, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.”

Table 1

declaration - electromagnetic emission	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies

Table 2

declaration - electromagnetic immunity		
Immunity test	IEC 60601 test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines	± 2 kV for power supply lines
Surge IEC 61000-4-5	± 0.5 kV, ± 1 kV line(s) to lines	± 0.5 kV, ± 1 kV line(s) to lines
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle At 0° , 45° , 90° , 135° , 180° , 225° , 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycles	0 % UT; 0.5 cycle At 0° , 45° , 90° , 135° , 180° , 225° , 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycles
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m
NOTE: UT is the a.c. mains voltage prior to application of the test level.		

Table 3

declaration - electromagnetic immunity		
Immunity test	IEC 60601 test level	Compliance level
Conducted RF IEC 61000-4-6	3 V 0.15 MHz to 80 MHz 6 V in ISM and amateur bands between 0.15 MHz and 80 MHz	3 V 0.15 MHz to 80 MHz 6 V in ISM and amateur bands between 0.15 MHz and 80 MHz
Radiated RF IEC 61000-4-3	10V/m 80 MHz to 2.7 GHz	10V/m

Table 4

declaration - IMMUNITY to proximity fields from RF wireless communications equipment					
Immunity test	IEC60601 test level				Compliance level
	Test frequency	Modulation	Maximum power	Immunity level	
Radiated RF IEC 61000-4-3	385 MHz	**Pulse Modulation: 18Hz	1.8W	27 V/m	27 V/m
	450 MHz	*FM+ 5Hz deviation: 1kHz sine	2 W	28 V/m	28 V/m
	710 MHz 745 MHz 780 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m
	810 MHz 870 MHz 930 MHz	**Pulse Modulation: 18Hz	2 W	28 V/m	28 V/m
	1720 MHz 1845 MHz 1970 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	2450 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	5240 MHz 5500 MHz 5785 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m
Note* - As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case. Note** - The carrier shall be modulated using a 50 % duty cycle square wave signal.					

Troubleshooting

Check the following if your unit should fail during operation. You can also refer to the pages of this manual for complete instructions.

Problem	Cause	Remedy
Nothing happens when the power switch is pressed.	Is the AC adapter connected correctly to the an electrical outlet and the compressor?	Check that the plug is inserted in an electrical outlet. Unplug then reinsert the plug if necessary.
No nebulisation or low nebulisation rate,when the power is on.	Is there medication in the medication tank?	Add the correct amount of medication to the medication tank.
	Is there too much/little medication in the medication tank?	
	Is the vaporiser head missing and not assembled correctly?	Attach the vaporiser head correctly.
	Is the nebulizer kit assembled correctly?	Assemble the nebulizer kit correctly.
	Is the nozzle blocked?	Make sure that the nozzle is free of blockages.
	Is the nebulizer kit tilted at a sharp angle?	Make sure that the nebuliser kit is not tilted at an angle of more than 45 degrees.
	Is the air tube connected correctly?	Make sure that the air tube is correctly connected to the compressor and nebuliser kit.
	Is the air tube folded or damaged?	Make sure that the air tube does not contain kinks.
	Is the air tube blocked?	Make sure that the air tube is free of blockages.
	Is the air filter dirty?	Replace the air filterwith a new one.

The device is abnormally loud.	Is the air filter cover attached correctly?	Attach the air filter cover correctly.
The device is very hot.	Is the compressor covered?	Do not cover the compressor with any type of cover during use.
	Is the device operating continuously over 20 minutes?	Limit use to 20 minutes at a time, and allow a 40 minutes interval before using the device again.

Note: If the suggested remedy does not solve the problem, do not try to repair the device- no parts of the unit are user serviceable.

Return the unit to an authorized retail outlet or distributor.

CONDITIONS FOR SERVICE UNDER GUARANTEE

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 2 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: Medicine Cup, Inhalation Masks, Mouthpiece, Nosepiece, Air Tube, Air Filters, DC Adapter
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@sanicare.it.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;
- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect use, installation or transport of the product;

- components or accessories subject to rapid wear and tear (electrode cups, elastic bands);
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and marketing purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

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Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

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Via/Piazza/ Street/Square: _____ N°: _____

Città/ City and State: _____ CAP/ Postal Code: _____

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DISPOSAL OF YOUR OLD APPLIANCE

When this crossed-out wheeled bin symbol is attached to a product, it means the product is covered by the European Directive 2002/96/EC.

All electrical and electronic products should be disposed off separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.

The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

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Shenzhen Bi-rich Medical Devices Co.,Ltd

Address:The 1st building of No. 10, Xinqiao GangZai Road, Xinqiao Street, Bao'An District ,518125, Shenzhen City, Guangdong Province, P. R. China

EC

REP

SUNGO Europe B.V.,

Add: Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel, The Netherlands

IMPORTATO DA

SANICO Srl

Via G.Ferraris, 37/41 - 20047 Cusago (MI) ITALY

info@sanicare.it - www.sanicare.it

Tel. (+39) 02 90 39 00 38

Fax (+39) 02 90 39 02 79

MADE IN CHINA

CE 1639