SANICO

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CONTENTS

1. Introduction and Intended Use
2. Important Information on Blood Pressure and its Measurement
3. Components of your blood pressure monitor
4. Using your Monitor for the First Time
5. Measurement Procedure
6. Care and Maintenance
7. Certification
8. Technical Specifications
9. EMC Declaration
10. Warranty

1.INTRODUCTION AND INTENDED USE

It enables reliable measurement of systolic and diastolic blood pressure as well as pulse through the oscillometric method.

Before using, please read this instruction manual carefully and then keep it in a safe place.

1.1. Remember...

- Only a health-care professional is qualified to interpret blood pressure measurements.
- This device is NOT intended to replace regular medical checkups.
- Blood pressure readings obtained by this device should be verified before prescribing or making adjustments to any medications used to control hypertension. Under no circumstances should YOU alter the dosages of any drugs prescribed by your physician.
- This monitor is intended for use by adults only. Consult with a physician before using this instrument on a child.
- In cases of irregular heartbeat, measurements made with this instrument should only be evaluated after consultation with a physician.
- Host products, including accessories, shall be processed in accordance with local regulations after reaching the life cycle.

1.2 Warnings and Precautions

Warning: The use of other accessories other than those specified or provided by the equipment manufacturer may cause electromagnetic radiation to increase or decrease electromagnetic immunity resulting in operational failure

Warning: This system may fail to yield specified measurement accuracy if operated or stored in temperature or humidity conditions outside the limits stated in the specifications section of this manual.

Warning: The separate ac adapter which is intended to connect USB interface of Blood Pressure Monitor has not been evaluated according to IEC 60601-1. The safety of the product shall be reappraised when it power supply by a separate ac adapter.

Warning: The user must check that the equipment functions safely and see that it is in proper working condition before being used.

Warning: No modification of this equipment is allowed.

Warning: The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.

Warning: The patient is an intended operator, the functions of monitoring blood pressure and pulse rate can be safely used by patient. The routine clean and changing batteries can be performed by the patient.

Warning: This device can not be used together with hf surgical equipment.

Warning: Use of power adapters

1. Adapter: input 100-240V, 50/60hz output d.c. 5V 1A
2. Do not to position the device to make it difficult to operate the disconnection device while using adaptor.
3. Do not be prone to water leakage, high temperature, moisture, direct sunlight and more or more corrosive gas environment. And Do not use this product in the above environment.

Warning: Too frequent measurements can cause injury to the PATIENT due to blood flow interference.

Warning: Don't place the cuff over wound part.

Warning: Pressurization of the CUFF can temporarily cause loss of function of simultaneously used monitoring ME EQUIPMENT on the same limb.

Caution: To avoid any possibility of accidental strangulation, keep this unit away from children and do not drape tubing around your neck.

Caution: To avoid damaging the device, keep this unit away from children and pets.

Caution: The standard material used for the bladder and tubing is latex-free.

Attention: Self-measurement means control, not diagnosis or treatment. Unusual values must always be discussed with your doctor. Under no circumstances should

you alter the dosages of any drugs prescribed by your doctor.

Attention: The pulse display is not suitable for checking the frequency of heart pacemakers!

Attention: In cases of irregular heartbeat, measurements made with this instrument should only be evaluated after consultation with a physician.

Note: To obtain the greatest accuracy from your blood pressure instrument, it is recommended that the instrument be used within the specified temperature and the relative humidity, please see the Technical Specifications

Note: The cuff is treated as the applied part. The user should contact the manufacturer for assistance, if needed, in setting up, using or maintaining the device.

Note: This device contains sensitive electronic components. Avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g. mobile telephones, microwave ovens) during use. These can lead to erratic results.

Note: Do not attempt to service or repair this device yourself. Should a malfunction occur, refer to local distributor or the manufacturer.

2. IMPORTANT INFORMATION ON BLOOD PRESSURE AND ITS MEASUREMENT

2.1. How does high or low blood pressure arise?


Your level of blood pressure is determined in the circulatory center of the brain and adjusts to a variety of situations through feedback from the nervous system. To adjust blood pressure, the strength and speed of the heart (Pulse), as well as the width of circulatory blood vessels is altered. Blood vessel width is controlled by fine muscles in the blood vessel walls.

Your level of arterial blood pressure changes periodically during heart activity: During the "blood ejection" (Systole) the value is highest (systolic blood pressure value). At the end of the heart's "rest period" (Diastole) pressure is lowest (diastolic blood pressure value).

Blood pressure values must lie within certain normal ranges in order to prevent particular diseases.


2.2. Which values are normal?

Please refer to the diagram below (Picture-01)

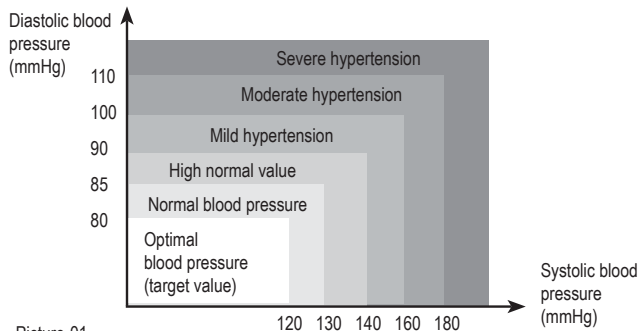Picture-01

There are six grids in the display of device. Please refer to the picture-01-01. Different grids represent different interval scales of WHO.

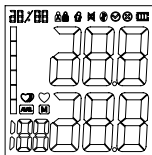
Blood pressure value	WHO grids in device	WHO Classification
DIA<80 & SYS <120	1	Optimal blood pressure
DIA<85 & SYS <130	2	Normal blood pressure
DIA<90 & SYS <140	3	High normal value
DIA<100 & SYS <160	4	Mild hypertension
DIA<110 & SYS <180	5	Moderate hypertension
DIA>=110 & SYS >=180	6	Severe hypertension


Picture-01-01

3. COMPONENTS OF YOUR BLOOD PRESSURE MONITOR

3.1. Measuring unit

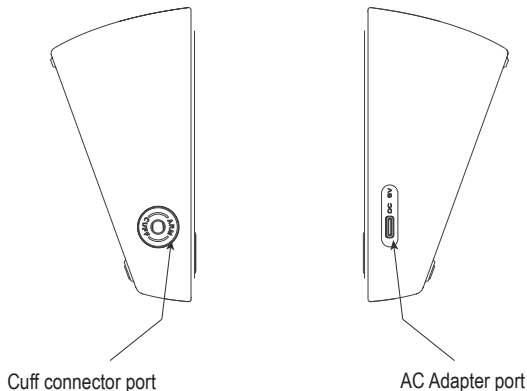Picture-02

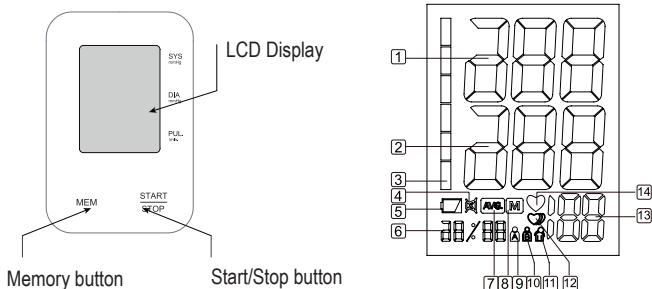Picture-03

3.2 The symbols on the LCD display

- | | |
|----------------------------|---|
| 1. Systolic blood pressure | 2. Diastolic blood pressure |
| 3. WHO function symbol | 4. Mute symbol |
| 5. Low battery symbol | 6. Date/Time display |
| 7. Average value symbol | 8. Memory symbol |
| 9. User A | 10. User B |
| 11. Movement error symbol | 12. Irregular heartbeat symbol |
| 13. Pulse symbol | 14. Heartbeat symbol (Flashes during measurement) |

3.3 Features of Model C03A


- | | |
|-----------------------------------|--------------------------------------|
| 1. Average value function | 2. Double users: 2 x 120 sets memory |
| 3. WHO function | 4. Irregular heartbeat checking |
| 5. External power adapter support | 6. Low battery display |
| 7. Date/time display | 8. Auto power-off |

4. USING YOUR MONITOR FOR THE FIRST TIME

4.1 Battery Installation

Use only 1.5V "AA" alkaline batteries with this device.

- Press the hook on the bottom of the battery cover and lift the cover off in the direction of the arrow.
- Install 3 "AA" size batteries so the + (positive) and (negative) polarities match the polarities of the battery compartment, then close the battery cover. Make sure that the battery cover is securely in position

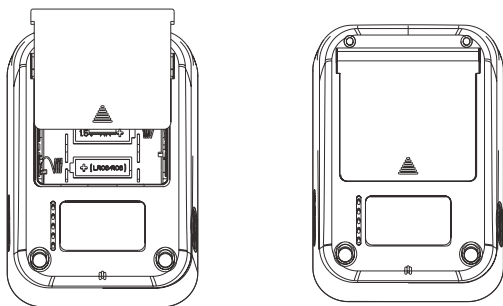Picture-04

Battery replacement

- Low Battery Indicator
- When the Low Battery Indicator appears on the display, turn the monitor off and remove all the batteries. Replace with 3 new batteries at the same time. Long-life alkaline batteries are recommended.
- To prevent the damage of monitor from leaked battery fluid, please take out of battery if the monitor unused in a long time (generally more than 3 months). If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Contact a physician immediately.
- Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution

4.2 System Settings

After you load the battery or connect power for the monitor, long press the START/STOP button for more than 3s, and then you can start to set.

Setting the User:

Press the MEM button to select User A or User B. When display A (/B) on the screen, press the MEM button to switch to user B (/A). Press the START/STOP button to confirm.

Setting the Year:

After setting the user, when the year display is flashing, press the MEM button continuously and it will increase continuously 1 by 1 until 2099, and then return the original year, once the year set is OK, press START/STOP button to confirm.

Setting Month/Date:

After setting the year, when the Month display is flashing, press the MEM button, the month will increase by 1, press START/STOP button to confirm, and do in the same way to set the date. Press START/STOP button to confirm.

Setting Time :

After setting the date, when the hour display is flashing, press the MEM button, the hour will increase by 1, press START/STOP button to confirm, and do in the same way to set the minute. Press START/STOP button to confirm.

Record Delete:

When you checking the memory data, long press MEM button to delete existing user measurement data.

Note:

please don't delete all the measurement records, if you decide to delete all the record, please keep the record in another way, in case you need it some days later. Take the battery out won't lead to a record missing.

4.2. Cuff tube connection

Insert the cuff tube into the opening on the left side of the monitor indicated by the drawing of a cuff.

5. MEASUREMENT PROCEDURE

5.1. Before measurement:

- Avoid eating and smoking as well as all forms of exertion directly before


measurement. These factors influence the measurement result. Find time to relax by sitting in an armchair in a quiet atmosphere for about ten minutes before taking a measurement.


- Remove any garment that fits closely to your upper arm.
- Always measure on the same arm (normally left).

5.3. Fitting the Cuff

Please refer to picture-05

- Wrap the cuff around your upper left arm. The rubber tube should be on the inside of your arm extending downward to your hand. Make certain the cuff lies approximately 2 to 3 cm above the elbow. Important! The **Artery Mark** on the edge of the cuff must lie over the artery which runs down the inner side of the arm.
- To secure the cuff, wrap it around your arm and press the hook and loop closure together.
- There should be little free space between your arm and the cuff. You should be able to fit 2 fingers between your arm and the cuff. Cuffs that don't fit properly result in false measurement values. Measure your arm circumference if you are not sure of proper fit.
- Lay your arm on a table (palm upward) so the cuff is at the same height as your heart. Make sure the tube is not kinked.

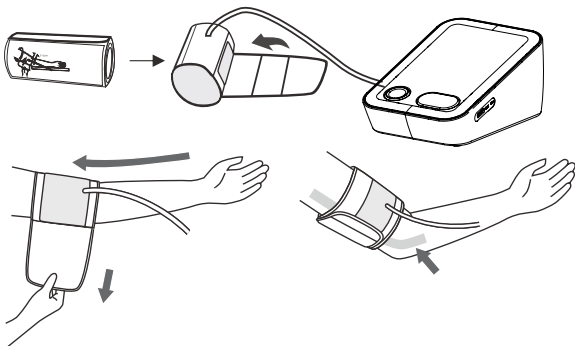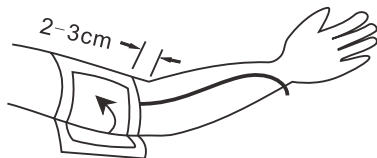Picture-05

5.4 Measure Procedure

Refer to picture 06

1. Sit comfortably in a chair with your feet flat on the floor.
2. Select the user ID (A or B): with the device off, press the START/STOP button for three seconds and change the user from A to B and vice versa, by pressing the MEM button; wait for the user symbol A or B to disappear from the display before proceeding with the measurement.
3. Stretch your arm forward on the desk and keep relaxing, make sure the palm of hand is upturned. Make sure arm is in correct position, to avoid body movement. Sit still and do not talk or move during the measurement.

After the cuff has been appropriately positioned on the arm and connected to the blood pressure monitor, the measurement can begin:

- a) Press the START/STOP button. The pump begins to inflate the cuff. In the display, the increasing cuff pressure is continually displayed.
- b) After automatically reaching an individual pressure, the pump stops and the pressure slowly falls. The cuff pressure is displayed during the measurement.
- c) When the device has detected your pulse, the heart symbol in the display begins to blink.
- d) When the measurement has been concluded, the measured systolic and diastolic blood pressure values, as well as the pulse will be displayed.
- e) The measurement results are displayed until you switch the device off if no button is pressed for 60 seconds, the device switches off automatically.

WARNING: THIS DEVICE MEASURES PRESSURE AND FREQUENCY DATA DURING INFLATION


- f) Movement error symbol ()

The Movement Error Symbol () is displayed if you move your body during the measurement. Please remove the cuff, and wait 2-3 minutes. Reapply the cuff and take another measurement.

NOTE:

Patient Position:

comfortably seated / legs uncrossed/ feet flat on the floor/ back and arm supported
middle of the CUFF at the level of the right atrium of the heart

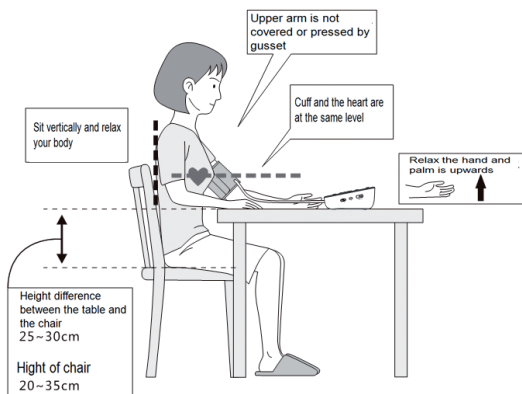Irregular Heartbeat Detector

This symbol  indicates that certain pulse irregularities were detected during the measurement.

In this case, the result may deviate from your normal basal blood pressure – repeat the measurement.

Information for the physician on frequent appearance of the Irregular Heartbeat Symbol

This instrument is an oscillometric blood pressure monitor device that also analyzes pulse frequency during measurement. The instrument is clinically tested.

If pulse irregularities occur during measurement, the irregular heartbeat symbol is displayed after the measurement. If the symbol appears more frequently (e.g. several times per week on measurements performed daily) or if it suddenly appears more often than usual, we recommend the patient to seek medical advice. The instrument does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

5.5. Error Indicates

SYMBOL	CAUSE	CORRECTION
No display appears	Weak battery or improper placement	Replace both batteries with new ones. Check the battery installation for proper placement of the battery polarities.
Er1	Sensor abnormal	Check if the pump is working or not. If it is working, then the problem is sensor abnormal. Please send it to the local distributor.
Er2	Monitor could not detect pulse wave or cannot calculate the blood pressure data	Check if the air releasing is too slow or not. If it is too slow, please check if there is any dust in the tube plug of the cuff and the cuff port in the device. If yes, please clean and start the measurement again. If no, please send the device back to the local distributor.
Er3	Measurement result is abnormal (SYS $\leq$ 35mmHg, DIA $\leq$ 23mmHg)	Occasionally-measure for one more time/ Always - send it to local distributor
Er4	Too loose cuff or air leakage (Cannot inflate to 30mmHg within 15s)	Tie the cuff correctly and make sure the air plug is properly inserted in the unit
Er 5	The air tube is crimped	Correct it and make the measurement again
Er 6	The sensor is sensing great fluctuation in the pressure	Please keep quiet and don't move
Er 7	The pressure that the sensor sensing is over the limit	Please send back to the local distributor
Er 8	The demarcation is incorrect or the device has not been demarcated	Please send back to the local distributor

Trouble removal

Problem	Check	Cause and solutions
No power	Check the battery power	Replace new ones
	Check the polarity position	Installation for proper placement of the batteries polarities
No Inflations	Whether the plug insert	Insert into the air socket tightly
	Whether the plug broken or leak	Change a new cuff
Err and stop working	Whether move the arm when inflate	Keep the body peaceful
	Check if chatting when measured	Keep quit when measure
Cuff leak	Whether the cuff wrap too lose	wrap the cuff tightly
	Whether the cuff broken	Change a new cuff

Please conotact the distributor if you can't solve the problem, do not disassemble the unit by yourself!

SYMBOL DESCRIPTIONS

The following symbols may appear in this manual, on the Digital Blood Pressure Monitor, or on it's accessories. Some of the symbols represent standards and compliances associated with the Digital Blood Pressure Monitor and its use.

Authorized Representative in the European Community

0123

CE Mark: conforms to essential requirements of the Medical Device Directive 93/42/EEC.

Date of manufacture

Manufacturer

Specifies serial number

Type BF applied part

Direct current

DISPOSAL: Do not dispose this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary.

Follow instructions for use

Medical Device

Put Up

Fragile

Afraid of the rain

Fear of the sun

Handle gently

Temperature range

No Sterilize requirement

Not category AP / APG equipment

Mode of operation: continuous

5.6. Memory

Each unit stores 120 sets measurements for 2 users, totally 240 sets (User A and B).

Viewing the stored values

With the unit off, press the Memory button. The display first shows "A", then shows an average of the 3 last measurements stored in the unit. Please note: Measurements for each user are averaged and stored separately. Be certain that you are viewing the measurements for the correct user. Pressing the Memory button again displays the previous value. To view a particular stored memory, press and hold the Memory button to scroll to that stored reading.

5.7. Discontinuing a Measurement

If it is necessary to interrupt a blood pressure measurement for any reason (e.g the patient feels unwell), the Power button can be pressed at any time. The device then immediately lowers the cuff pressure automatically.

5.10. Using the AC Adapter

You may also operate this monitor using the AC adapter (output 5V DC/1A with Type-C connector).

Use only the approved AC adapter to avoid damaging the unit.

- a) Ensure that the AC adapter and cable are not damaged.
- b) Plug the adapter cable into the AC adapter port on the right side of the blood pressure monitor.
- c) Plug the adapter into your electrical outlet. When the AC adapter is connected, no battery current is consumed.

Note: No power is taken from the batteries while the AC adapter is connected to the monitor. If electrical power is interrupted, (e.g., by accidental removal of the AC adapter from the outlet) the monitor must be reset by removing the plug from the socket and reinserting the AC adapter connection.

6. CARE AND MAINTENANCE

Wash hands after each time measurement.

If one device is used by different patients, wash hands before and after each use.

- a) Do not expose the device to either extreme temperatures, humidity, dust or direct sunlight.
- b) The cuff contains a sensitive air-tight bubble. Handle this cuff carefully and avoid all types of stress through twisting or buckling.
- c) Clean the device with a soft, dry cloth. Do not use gas, thinners or similar solvents. Spots on the cuff can be removed carefully with a damp cloth and soapsuds. The cuff with bladder must not be washed in a dishwasher, clothes washer, or submerged in water.
- d) Handle the tube carefully. Do not pull on it. Do not allow the tubing to kink and keep it away from sharp edges.
- e) Do not drop the monitor or treat it roughly in any way. Avoid strong vibrations.
- f) Never open the monitor! This invalidates the manufacturer's warranty.
- g) Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

6.1. Accuracy test

Sensitive measuring devices must be checked for accuracy from time to time. We recommend a periodical inspection of your unit by an authorized dealer every 1 years. Please turn to local distributor or the manufacturer.

7. CERTIFICATIONS

Device standard:

This device is manufactured to meet the European blood pressure monitors:

EN1060-3 / IEC 80601-2-30 / ISO81060-1 / IEC60601-1-11 / IEC60601-1

Electromagnetic compatibility:

Device fulfills the stipulations of the International standard

IEC60601-1-2

8. TECHNICAL SPECIFICATIONS

Model: PL094(C03A)

Weight: 164g (weight of the device, batteries and AC adapter not included)

Display: 66mm x49,2mm LCD Digital Display

Size: 123 (W) x 85 (L) x 64(H) mm

Accessories: 1×Main Device, 1×Cuff, 1×Users manual

Operating Conditions: Temperature: 5° to 40°; Humidity: 15% to 93% RH;

Pressure altitude: 70KPa~ 106Kpa

Storage And Shipping Conditions: Temperature: -25° to 70°;

Humidity: ≤ 93% RH; Atmospheric Pressure range: 70KPa~ 106Kpa

Measuring method: Oscillometric

Pressure sensor: Restive

Measuring range: DIA:40-130mmHg; SYS: 60-230mmHg

Pulse: 40 to 170 per minute

Cuff pressure display range: <300 mmHg

Memory: Automatically stores the last 120 measurements for 2 users (total 240)

Measuring resolution: 1 mmHg

Accuracy: Pressure within ± 3 mmHg / pulse ± 5 % of the reading

Power source: a) 3* AA batteries, 1.5 V

b) AC adapter a.c. 100-240V 50/60HZ OUTPUT: d.c. 5V 1A

Accessories: Wide range rigid cuff : 8.7" – 15.7" (22 - 40 cm)

Automatically power off : 60 seconds

Users: Adult

Expected service life of the device and accessories: 5 years

Technical alterations reserved!

9. EMC Declaration

- 1) *This product needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided, and this unit can be affected by portable and mobile RF communications equipment.
- 2) * Do not use a mobile phone or other devices that emit electromagnetic fields, near the unit. This may result in incorrect operation of the unit.
- 3) * Caution: This unit has been thoroughly tested and inspected to assure proper performance and operation!
- 4) * Caution: this machine should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, this machine should be observed to verify normal operation in the configuration in which it will be used

Guidance and manufacturer's declaration – electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should ensure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrical fast transient/burst IEC 61000-4-4	+/-2 kV for power supply lines +/-1 kV for input/output lines	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	+/-1 kV differential mode +/-2 kV common mode	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0 % UT ; 1 cycle 70 % UT; 25/30 cycle 0% UT; 250/300 cycle	Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m 50/60HZ	3 A/m 50/60HZ	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE UT is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration – electromagnetic immunity

The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should ensure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz 3 V RMS outside the ISM band, 6 V RMS in the ISM and amateur bands 80% AM at 1kHz	Not applicable	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 0.35\sqrt{P}$ $d = 1.2\sqrt{P}$
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz 80% AM at 1kHz	10 V/m 80 MHz to 2.7 GHz 80% AM at 1kHz	80MHz to 800MHz: $d = 1.2\sqrt{P}$ 800MHz to 2.7GHz: $d = 2.3\sqrt{P}$ Where, P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance. Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the device.

Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Guidance and manufacturer's declaration – electromagnetic emission		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should ensure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The "Digital Blood Pressure Monitor PL096" uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Not applicable	

Recommended separation distances between portable and mobile RF communications equipment and the device			
The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1,2\sqrt{P}$	80 MHz to 800 MHz $d = 1,2\sqrt{P}$	800 MHz to 2,5 GHz $d = 2,3\sqrt{P}$
0,01	0.12	0.12	0.23
0,1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Guidance and manufacturer's declaration - electromagnetic immunity						
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device, should assure that it is used in such an environment.						
Test frequency (MHz)	Band a) (MHz)	Service a)	Modulation a)	Maximum power (w)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
385	380-390	TETRA 400	Pulse Modulation b) 18 Hz	1.8	0.3	27
450	430-470	GMRS 460, FRS 460	FM c) ± 5 kHz deviation 1 kHz sine	2	0.3	28
710	704-787	LTE Band 13, 17	Pulse Modulation b) 217 Hz	0.2	0.3	9
745						
780						
810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse Modulation b) 18 Hz	2	0.3	28
870						
930						
1720						
1845	1700-1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1,3 4,25;UMTS	Pulse Modulation b) 217 Hz	2	0.3	28
1970						
2450	2400-2570	Bluetooth, WLAN 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse Modulation b) 217 Hz	2	0.3	28
5240	5100-5800	WLAN 802.11 a/n	Pulse Modulation b) 217 Hz	0.2	0.3	9
5500						
5785						

NOTE If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

a) For some services, only the uplink frequencies are included.

b) The carrier shall be modulated using a 50% duty cycle square wave signal.

c) As an alternative to FM modulation. 50% pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

The MANUFACTURER should consider reducing the minimum separation distance, based on RISK MANAGEMENT, and using higher IMMUNITY TEST LEVELS that are appropriate for the reduced minimum separation distance. Minimum separation distances for higher IMMUNITY TEST LEVELS shall be calculated using the following equation:

$$E = \frac{6}{d} \sqrt{P}$$

Where P is the maximum power in W, d is the minimum separation distance in m, and E is the IMMUNITY TEST LEVEL in V/m.

9. CONDITIONS FOR SERVICE UNDER GUARANTEE

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 2 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: electrode cups, elastic bands, power supply cables etc..
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@sanicare.it.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;
- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect use, installation or transport of the product;

- components or accessories subject to rapid wear and tear (electrode cups, elastic bands);
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and profit-making purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

e-mail: info@sanicare.it

www.sanicare.it

Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

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SECURITY OF ALKALINE BATTERIES INTO PRODUCTS

1. Fire - Never delete the batteries throwing them in a fire. they could break and release the substances contained inside them.
2. If ingested or inserted, seek immediate medical assistance.
3. Do not store or leave the batteries in places accessible by children.
4. Remove old or weak batteries from appliances.
5. Do not mix old and new batteries in the same apparatus. There may be losses. The replacement of all batteries in a device must be done simultaneously.
6. Never mix different types of batteries in a single device. There may be losses.
7. Install carefully the batteries in the correct polarity (+ and -) according to the instructions provided by the manufacturer of the product
8. Do not leave the Products turn on with batteries when they are depleted. There may be losses.
9. Remove the batteries from any device that must remain unused for long periods of time
10. Avoid contact of the batteries with metal objects or other batteries because they may experience losses or they could break.
11. Do not recharge batteries because there may be some losses or they could break.
12. Not groped to disassemble the batteries
13. Store batteries in a cool and dry. Avoid sudden changes in temperature. Store the batteries in their original packaging before use.
14. RESPECT THE ENVIRONMENT - DISPOSE OF EXHAUSTED BATTERIES IN THE CONTAINERS

The battery and the possible Electronic Equipment -AEE- or anything that contains it must not be disposed as urban waste, but in compliance with Directive 2002/96/EC, 2006/66/EC and of National Laws transposition for the disposal of the product end of life, and must be subjected to a separate collection to avoid damage to the environment. Failure to comply with the above is sanctioned by Law.

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DISPOSAL OF YOUR OLD APPLIANCE

When this crossed-out wheeled bin symbol is attached to a product, it means the product is covered by the European Directive 2002/96/EC.

All electrical and electronic products should be disposed off separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.

The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

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WARRANTY IS VALID ONLY IF ACCOMPANIED BY INVOICE/TICKET.

Company name: Shenzhen Jamr Technology Co., Ltd.

Address: A101-301, D101-201, Jamr Science & Technology Park, No. 2
Guiyuan Road, Guixiang Community, Guanlan Street, Longhua District,
Shenzhen 518100, PEOPLE'S REPUBLIC OF CHINA

Company Name: Shanghai International Holding Corp. GmbH (Europe)

Address: Eiffestrasse 80, 20537 Hamburg, Germany
shholding@hotmail.com

IMPORTATO DA

SANICO Srl

Via G.Ferraris, 37/41 - 20047 Cusago (MI) ITALY

info@sanicare.it - www.sanicare.it

Tel. (+39) 02 90 39 00 38

Fax (+39) 02 90 39 02 79

Versione: 5022 / Mod. PL094 C03A

MADE IN CHINA

C € 0123