SANICO+

MISURATORE DI PRESSIONE
TENSIOMÈTRE
BLUTDRUCKMESSGERÄT
BLOOD PRESSURE MONITOR
BLOEDDRUKMETER
TENSIÓMETRO

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Signature:

Date: / /

WARRANTY IS VALID ONLY IF ACCOMPANIED BY INVOICE/TICKET.

Table of Contents

1. Introduction

2. Background information on blood pressure

- A. Illustration of the blood pressure monitor
- B. Appearance of the Irregular Heartbeat (IHB) Symbol

3. Putting the blood pressure monitor into operation

- A. Battery installation
- B. AC/DC Adapter

4. How to take blood pressure measurement

- A. Preparation before BP measurement
- B. Evaluation of my blood pressure
- C. Attaching the arm cuff
- D. Measurement procedure
- E. Displaying the last measurement
- F. Discontinuing a measurement

5. Troubleshooting guide

- A. Error warnings
- B. What to do if the unit does not work

6. Safety, Care, Accuracy Test

7. Specifications

8. Warranty

1. Introduction

Congratulations on your purchase of the PL097 blood pressure monitor. This unit uses a new oscillometric method to provide rapid measurement of systolic and diastolic pressure as well as pulse rate.

This unit has hospital-grade accuracy, and has a jumbo LCD display to clearly show the operating status of cuff pressure and pumping pressure. These unique features allow for easy recognition of pressure changes during measurement.

Before using this unit, be sure to read all the instructions very carefully. Please contact your physician or your pharmacist if you have any questions.

Intended use:


This oscillometric blood pressure monitor is intended for measuring non-invasive blood pressure in people aged 12 years or older.

It is clinically validated in patients with hypertension, hypotension, diabetes, pregnancy, pre-eclampsia, atherosclerosis, end-stage renal disease, obesity and the elderly.

2. Background information on blood pressure

A. Illustration of the blood pressure monitor

The illustration shows the Automatic Blood Pressure Monitor with "One Touch" operation consisting of:

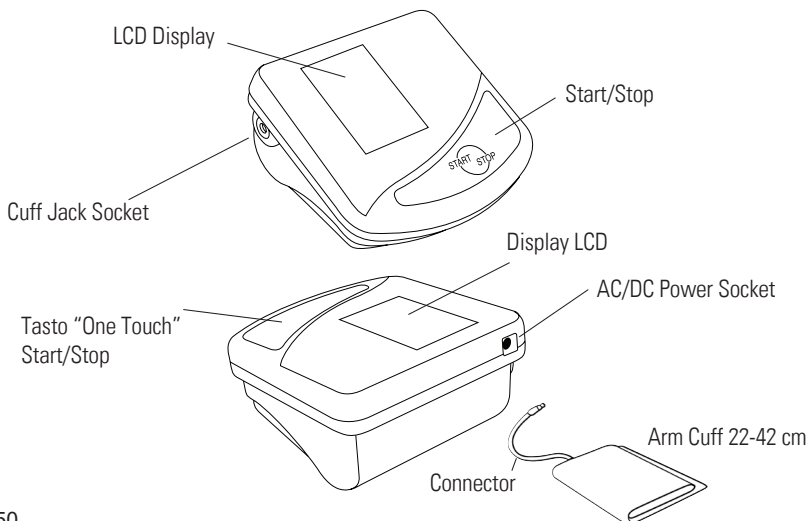B. Appearance of the Irregular Heartbeat (IHB) Symbol

This symbol  indicates that an irregular heartbeat was detected. In this case, the measured blood pressure may deviate from your actual blood pressure values. It is recommended to repeat the measurement.

Information for the doctor in case of repeated appearance of the IHB symbol:

This device is an oscillometric blood pressure monitor that also measures the pulse during blood pressure measurement and indicates when the heart rate is irregular.

3. Putting the blood pressure monitor into operation.

A) Inserting the batteries

After you have unpacked your device, first insert the batteries. The battery compartment is located on the underside of the device (see illustration).


a) Remove cover as illustrated.

b) Insert batteries (Use 4 new, long-life 1.5 V, size AA alkaline batteries), observing the indicated polarity.

c) If the battery warning appears in the display, the batteries are low and must be replaced by new ones.


d) Do not mix old and new batteries.

Do not mix alkaline, standard (carbon-zinc), or rechargeable (ni-cad, ni-mh, etc) batteries.

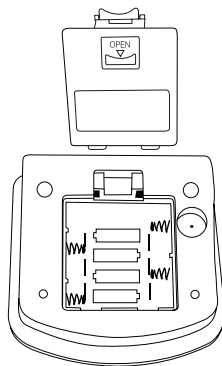Attention!

- After the battery warning appears, the device is blocked until the batteries have been replaced.
- Use 4 new, long-life 1.5 V, size AA alkaline batteries.
- Do not use batteries beyond their date of expiry.
- Remove batteries if the device is not going to be used for a prolonged period.

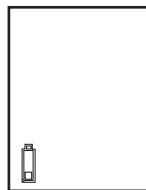Functional check: Hold the Start/Stop button down to test all the display elements. When functioning correctly all segments must appear.

Insert the connector on the tube cuff into the cuff jack provided on the left side of the instrument

Tube connection

(Tube cuff + connector):

Set up (or check) the tube connections between cuff and the instrument.

B) Using an AC/DC adaptor (special accessory - not included)

It is possible to operate this blood pressure instrument, through an AC/DC adaptor not included (output 6 V DC / 600mA with DIN plug). Make certain that you use an AC/DC adaptor which fulfils the legal requirements. These can be purchased at most electronics stores. Ensure that neither the mains adapter nor the cable are damaged.

The AC/DC adaptor should not exceed 600 mA


a) Plug the DIN plug into the socket at the back of the instrument.

b) Plug the AC/DC adaptor into a 120 V power socket.

Test that power is available by pressing the Start/Stop button.

Note:

- No power is taken from the batteries while the AC/DC adaptor is connected to the instrument.
- If the voltage is interrupted during the measurement (e.g. by removal of the AC/DC adaptor from the socket), the instrument must be reset by removing the plug from its socket.

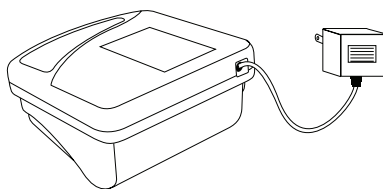4. How to take a measurement

A) Preparation

- Avoid activity, eating or smoking immediately before the measurement.
- Sit down on a back-supported chair and relax for 5 minutes.
- Keep the feet flat on the floor and do not cross your legs.
- Always measure on the same arm (normally left). It is recommended that doctors perform double arm measurements on a patients first visit in order to determine which arm to measure in the future. The arm with the higher blood pressure should be measured.
- Remove close-fitting garments from the upper arm. To avoid constriction, shirt sleeves should not be rolled up - they do not interfere with the cuff if they are laid flat.
- Always ensure that the correct cuff size is used (marking on the cuff).
- Fit the cuff closely, but not too tight.
- Make sure that the cuff is positioned 1-2 cm above the elbow.
- The artery mark on the cuff (ca. 3 cm long bar) must lie over the artery which runs down the inner side of the arm.
- Support your arm so it is relaxed.
- Ensure that the cuff is at the same height as your heart.

B) Evaluation of my blood pressure

The classification corresponds to the following ranges defined by international guidelines (ESH, ESC, JSH). Data in mmHg.

Range		Systolic	Diastolic	Recommendation
1.	blood pressure too high	≥135	≥85	Seek medical advice
2.	“blood pressure elevated”	130 - 134	80 - 84	Self-check
3.	blood pressure normal	<130	<80	Self-check

The higher value is the one that determines the evaluation.


Example: a blood pressure value of 140/80 mmHg or a value of 130/90 mmHg indicates «blood pressure too high».

C) Attaching the arm cuff



1. Pass the end of the cuff through the flat metal ring so a loop is formed. The closure must be facing outward. (Ignore this step if the cuff has already been prepared).
2. Place the cuff around your arm. Make sure the bottom edge of the cuff is about 1 to 2 cm above the elbow joint. Adjust the cuff so that the rubber tubing under the cuff lies over the brachial artery, which runs on the inside of the arm (see Fig. C).
3. Tighten the cuff by pulling the end and close the cuff.
4. The cuff should fit snugly around the arm, but not too tight.
5. Place the arm on the table (palm facing upwards) so that the cuff is at the same level as the heart. Make sure there is no kink in the hose, and your arm is resting on the table.
6. You can adjust the level of your arm by putting a cushion under your arm.
7. Remain seated in a comfortable room temperature for at least 5 minutes, then start the measurement.
8. For those who cannot put the cuff on the left arm, put it on the right arm as shown.
9. Consecutive measurements will cause blood accumulation in the lower arm which will affect the measuring results. To improve reading accuracy, raise the arm being measured, squeeze and relax your hand several times, then take another measurement. Another option is to take the cuff off and wait a few minutes before repeating measurement.

Cuff
Assembly

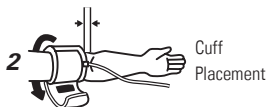Cuff
Placement

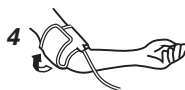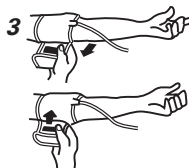
Cuff at
Heart Level

Adjust Arm
Level as
Needed

Sit Quietly
During
Measurement

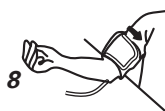D) Measurement Procedure

After the cuff has been appropriately positioned, the measurement can begin:

a) Press the Start/Stop button. The pump begins to inflate the cuff. The rising pressure in the cuff is shown on the display.


b) After the suitable inflation pressure has been reached, the pump stops and the pressure gradually falls. The cuff pressures are displayed.

c) When the instrument detects a pulse, the heart symbol on the display starts to flash and a beep is heard for every heart beat.

d) A longer beep is sounded when the measurement has been completed. The systolic (top), diastolic (middle) blood pressures and pulse rate (bottom) now appear in the LCD display.

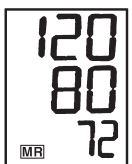
e) The results of the measurement remain on the display until you switch off the instrument. If no button is pressed for a period of 5 minutes the instrument switches itself off in order to preserve the batteries.

E) Memory – displaying the last measurement

The measured results are stored in the instrument until a new measurement is taken or the batteries are removed. You can call up the stored values at any time by pressing the Start/Stop button for at least 5 seconds.

F) Discontinuing a measurement

If it is necessary to interrupt a blood pressure measurement for any reason (e.g. the patient feels unwell), the Start/Stop button can simply be pressed at any time. The instrument then immediately releases the cuff pressure.

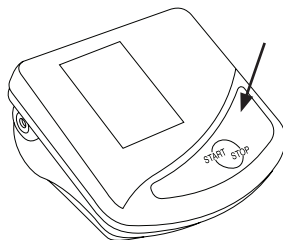5. Troubleshooting guide

A) Error warnings

Error warnings will appear on the LCD, as shown, if the following occur:

Error	Possible reasons and solutions
Err 1	The measurement of the pressure is impossible. The tube of the cuff is not connected correctly or is folded. Please verify the connection and the tube of the cuff.
Err 2	The inflation is insufficient or the user has moved or talked during the measurement. Please stay calm for one minute and restart the measurement without moving and talking.
Err 3	The inflation of the cuff is incorrect. Please verify that the cuff is properly wrapped and that the connection of the tube is correct. Restart the measurement.
Err 5	The difference between systolic and diastolic is excessive. Measure again carefully following proper cuff procedures and ensure measurement under quiet conditions.
HI	The pressure in the cuff is too high (over 300 mmHg) OR the pulse is too high (over 200 beats per minute). Relax for 5 minutes and repeat the measurement.*
LO	The pulse is too low (less than 40 beats per minute). Repeat the measurement.*

* Please consult your doctor, if this or any other problem occurs repeatedly.

B) How to troubleshoot (Unit does not work)

Trouble	Correction
Abnormal display or LCD blank when power is turned on, though batteries are installed.	1) Check and correct the polarity of the installed batteries. 2) Re-install the batteries or replace with new ones.
The unit does not measure the blood pressure values or abnormally low/high values are displayed.	1) Wrap the cuff in correct position again. 2) Sleeve may be rolled up and constricting arm. Remove clothing. 3) Measure blood pressure again.
Please Note: If a patient has serious irregular heartbeat, it may not always be possible to measure blood pressure.	
Measured values are different from those measured in hospital; or, every measurement is different.	1) Blood pressure always changes because values depend on physical mental conditions during the time of measurement. 2) Record a trend of daily change and consult with your physician for guidance. 3) Many patients experience elevated readings in their physician's office due to apprehension.
The pressure does not rise although the pump is running.	Check the connection of the cuff tube and connect properly if necessary.

Note: Some measuring errors result from special physical conditions and some may be related to user technique; but in rare cases, there can also be mechanical problems. If the unit still does not work, even after the above, bring it to your pharmacist. Do not disassemble the unit and attempt to repair it yourself. (There are no user serviceable parts inside).

ADDITIONAL INFORMATION:

Blood pressure is also subject to variations in healthy people. However, it must be stressed that measurements can only be compared if they were made under identical conditions (resting, and at same time of day)! If the variations are greater than 15 mmHg in spite of observing all these factors, and / or you hear an irregular pulse on several occasions, please consult your physician. You should consult your pharmacist if there are technical problems with the unit. Any unauthorized repairs will invalidate your guarantee.

6. Safety, Care, Accuracy Test

Safety and protection

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device may only be used for the purposes described in these instructions. The manufacturer cannot be held liable for damage caused by incorrect application.
- This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in the «Technical Specifications» section.
- Protect it from:
 - water and moisture /extreme temperatures/ impact and dropping/ contamination and dust/ direct sunlight /heat and cold
- The cuffs are sensitive and must be handled with care.
- Do not exchange or use any other kind of cuff or cuff connector for measuring with this device.
- Only pump up the cuff once fitted.
- Do not use this device close to strong electromagnetic fields such as mobile telephones or radio installations. Keep a minimum distance of 3.3 m from such devices when using this device.
- Do not use this device if you think it is damaged or notice anything unusual.
- Never open this device.
- If the device is not going to be used for a prolonged period the batteries should be removed.
- Read the additional safety information provided within the individual sections of this instruction manual.
- If you suffer from a cardiac arrhythmia consult with your doctor before using the device. See also chapter «Appearance of the Irregular Heartbeat (IHB) Symbol» of this user manual.
- The measurement results given by this device is not a diagnosis. It is not replacing the need for the consultation of a physician, especially if not matching the patient's symptoms. Do not rely on the measurement result only, always consider other potentially occurring symptoms and the patient's feedback. Calling a doctor or an ambulance is advised if needed.
- Permanently high blood pressure values can damage your health and must be treated by your doctor!
- Always discuss your values with your doctor and tell him/her if you have noticed anything unusual or feel unsure. Never rely on single blood pressure readings.
- Under no circumstances should you alter the dosages of drugs or initiate a treatment without consulting your doctor.
- Deviations between measurements taken by your doctor or in the pharmacy and those taken at home are quite normal, as these situations are completely different.
- The pulse display is not suitable for checking the frequency of heart pacemakers!
- If you are pregnant, you should monitor your blood pressure regularly as it can change drastically during this time.

Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed. Be aware of the risk of strangulation in case this device is supplied with cables or tubes.

Device care

- Clean the device only with a soft, dry cloth.

Cleaning the cuff

Carefully remove spots on the cuff with a damp cloth and soapsuds.

WARNING: Do not wash the cuff in a washing machine or dishwasher!

Accuracy test

- We recommend this device is tested for accuracy every 2 years or after mechanical impact (e.g. being dropped). Please contact your local Service to arrange the test.

7. Specifications

Weight:	410 g (with batteries)
Size:	146 (W) x 145 (L) x 46 (H) mm
Storage Temperature:	-20°C to 55°C (-4°F to 131°F)
Storage Humidity:	15 to 95% relative humidity maximum
Operating Temperature:	10°C to 40°C (50°F to 104°F)
Operating Humidity:	15 to 95% relative humidity maximum
Indicator:	LCD (Liquid Crystal Display)
Measuring Method:	Oscillometric method corresponding to Korotkoff method: Phase I systolic, Phase V diastolic
Measuring range:	30 to 280 mmHg
SYS/DIA:	40 to 200 per minute
Pulse:	
Cuff pressure display range:	0 - 299 mmHg
Resolution:	1 mmHg
Memory:	Stores the last set of measurements automatically
Static accuracy:	Pressure within ± 3 mmHg
Pulse accuracy:	± 5 % of the readout value
Power source:	a) 4 x 1.5 V alkaline batteries, size AA b) mains adapter 6 V DC 600mA (voltage 4.5 V DC to 6 V DC)
Cuff:	arm cuff circumference 22-42 cm
Device standard:	Device corresponds to the requirements of the standard for non-invasive blood pressure monitor EN 1060-1 EN 1060-3 EN 1060-4 IEC 60601-1 IEC 60601-1-2 IEC 60601-1-11
Electromagnetic compatibility:	Device fulfils the stipulations of the standard IEC 60601-1-2

Battery lifetime: approx. 800 measurements (using new batteries)

Expected service life: Device: 5 years or 10000 measurements

Accessories: 2 years

This device complies with the requirements of the Medical Device Directive 93/42/EEC.

Technical alterations reserved.

8. CONDITIONS FOR SERVICE UNDER GUARANTEE

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 5 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: electrode cups, elastic bands, power supply cables etc. nor regular checks and maintenance (calibration)
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@sanicare.it.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;

- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect use, installation or transport of the product;
- components or accessories subject to rapid wear and tear (electrode cups, elastic bands);
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and profit-making purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

e-mail: info@sanicare.it

www.sanicare.it

Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

Send this coupon in case of repairs:

Model: PL097

Serial Number:

(if available)

Date of purchase: / /

Dealer's Stamp:

Buyer's Full Name:

Street/Square: N°:

City and State: Postal Code:

Country:

Phone Number:

E-mail: @

Problem description:

.....
.....
.....

Signature:

Date: / /

WARRANTY IS VALID ONLY IF ACCOMPANIED BY INVOICE/TICKET.

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DISPOSAL OF YOUR OLD APPLIANCE

When this crossed-out wheeled bin symbol is attached to a product, it means the product is covered by the European Directive 2002/96/EC.

All electrical and electronic products should be disposed off separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.

The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

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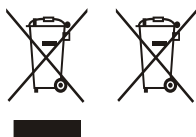
SECURITY OF ALKALINE BATTERIES INTO PRODUCTS

1. Fire - Never delete the batteries throwing them in a fire. they could break and release the substances contained inside them.
2. If ingested or inserted, seek immediate medical assistance.
3. Do not store or leave the batteries in places accessible by children.
4. Remove old or weak batteries from appliances.
5. Do not mix old and new batteries in the same apparatus. There may be losses. The replacement of all batteries in a device must be done simultaneously.
6. Never mix different types of batteries in a single device. There may be losses.
7. Install carefully the batteries in the correct polarity (+ and -) according to the instructions provided by the manufacturer of the product
8. Do not leave the Products turn on with batteries when they are depleted. There may be losses.
9. Remove the batteries from any device that must remain unused for long periods of time
10. Avoid contact of the batteries with metal objects or other batteries because they may experience losses or they could break.
11. Do not recharge batteries because there may be some losses or they could break.
12. Not groped to disassemble the batteries
13. Store batteries in a cool and dry. Avoid sudden changes in temperature. Store the batteries in their original packaging before use.
14. RESPECT THE ENVIRONMENT - DISPOSE OF EXHAUSTED BATTERIES IN THE CONTAINERS

The battery and the possible Electronic Equipment -AEE- or anything that contains it must not be disposed as urban waste, but in compliance with Directive 2002/96/EC, 2006/66/EC and of National Laws transposition for the disposal of the product end of life, and must be subjected to a separate collection to avoid damage to the environment. Failure to comply with the above is sanctioned by Law.

Microlife Corporation
9F, 431, Ruiguang Road, Neihu
Taipei 11492
Taiwan, China

Microlife UAB
P. Lukšio g. 32
08222 Vilnius, Lithuania

Microlife AG
Esenstrasse 139
9443 Widnau / Switzerland

MPORTATO DA
SANICO Srl
Via G.Ferraris, 37/41 - 20047 Cusago (MI) ITALY
info@sanicare.it - www.sanicare.it
Tel. (+39) 02 90 39 00 38
Fax (+39) 02 90 39 02 79

C € 0044

MADE IN CHINA