
SANICO+

MISURATORE DI PRESSIONE
BLOOD PRESSURE MONITOR
TENSIOMÈTRE
BLOEDDRUKMETER
BLUTDRUCKMESSGERÄT
TENSIÓMETRO

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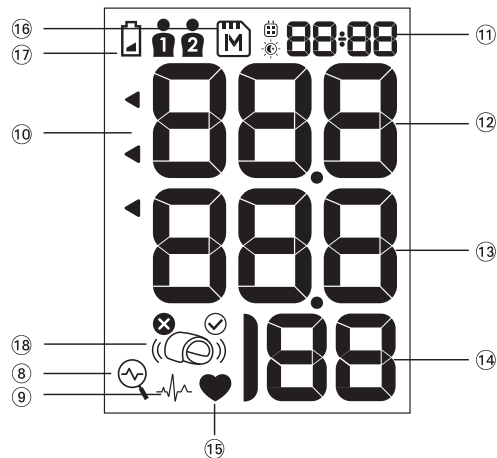- ① ON/OFF button
- ② Display
- ③ M-button (memory)
- ④ TIME-button
- ⑤ Traffic Light
- ⑥ Battery Compartment
- ⑦ Cuff

Display

- ⑧ Cuff Signal Indicator «Err 1»
- ⑨ Irregular Heartbeat Detection
- ⑩ Traffic Light Indicator
- ⑪ Date/Time
- ⑫ Systolic Value
- ⑬ Diastolic Value
- ⑭ Pulse Rate
- ⑮ Pulse Indicator
- ⑯ Stored Value
- ⑰ Battery Display
- ⑱ Cuff Fit Check
- A: Suboptimal Cuff Fit
- B: Arm Movement Indicator «Err 2»
- C: Cuff Pressure Check «Err 3»

Intended Use, User, and Environment:

This device is an oscillometric blood pressure monitor intended for radial (wrist) blood pressure non invasively.

It is suitable for use by adults (age 12 and above) for self-measurement and measurement of others, in a domestic environment (e.g. home-use).

Dear Customer,

This monitor enables very fast and reliable measurement of the systolic and diastolic blood pressure as well your pulse by use of oscillometric method. This device was developed in collaboration with physicians and clinical tests carried out prove its measurement accuracy to be of a very high standard*.

If you have any question, problems or want to order spare parts your local distributor's service. Your dealer or pharmacy will be able to give you the address of your local distributor's. Alternatively, visit the internet at www.sanicare.it where you will find a complete information on our products.

** This device is tested according to the ESH.*

Read the instructions carefully before using this device.

Medical Device - Wrist Blood Pressure Monitor

Type BF applied part

Keep Dry

Table of Contents

1. Important Facts about Blood Pressure and Self-Measurement

2. Using the Device for the First Time

- Inserting the batteries
- Setting the date and time

3. Checklist for Taking a Reliable Measurement

4. Taking a Blood Pressure Measurement

- How not to store a reading
- How do I evaluate my blood pressure?

5. Appearance of the Irregular Heart Beat (IHB)

6. Data Memory

- Viewing the stored values
- Clearing all values

7. Battery Indicator and Battery change

- Low battery
- Flat battery replacement
- Which batteries and which procedure?
- Using rechargeable batteries

8. Error Messages

9. Safety, Care, Accuracy Test and Disposal

- Safety and protection
- Device care
- Cleaning the cuff
- Accuracy test
- Disposal

10. Conditions for service under guarantee

11. Technical Specifications

Guarantee Card (see Back Cover)

1. Important Facts about Blood Pressure and Self-Measurement

Blood pressure is the pressure of the blood flowing in the arteries generated by the pumping of the heart. Two values, the systolic (upper) value and the diastolic (lower) value, are always measured.

The device indicates the pulse rate (the number of times the heart beats in a minute). Permanently high blood pressure values can damage your health and must be treated by your doctor!

Always discuss your values with your doctor and tell him/her if you have noticed anything unusual or feel unsure. Never rely on single blood pressure readings.

There are many causes of excessively high blood pressure values. Your doctor will explain them in more detail and offer treatment where appropriate. Besides medication, relaxation technique, weight loss and exercise can also lower your blood pressure.

Under no circumstances should you alter the dosages of any drugs prescribed by your doctor!

Depending on physical exertion and condition, blood pressure is subject to wide fluctuations as the day progresses. You should therefore take your measurements in the same quiet conditions and when you feel relaxed! Take at least two measurements per day, one in the morning and one in the evening.

It is quite normal for two measurements taken in quick succession to produce significantly different results.

Deviations between measurements taken by your doctor or in the pharmacy and those taken at home are quite normal, as these situations are completely different.

Several measurements provide a much clearer picture than just one single measurement.

Leave a small break of at least 15 seconds between two measurements.

If you are pregnant, you should monitor your blood pressure very closely as it can change drastically during this time!

If you suffer from an irregular heartbeat, measurements taken with this device should only be evaluated after consultation with your doctor.

The pulse display is not suitable for checking the frequency of heart pacemakers!

2. Using the Device for the First Time

Inserting the batteries

After you have unpacked your device, first insert the batteries. The battery compartment ⑥ is on the bottom of the device. Insert the batteries (2x1.5V, size AAA), thereby observing the indicated polarity.

Setting the date and time

1. After the new batteries are fitted, the year number flashes in the display. You can set the year by pressing the M-button ③ to confirm and then set the month, press the TIME-button ④
2. Press the M-button to set the month. Press the TIME-button to confirm and then set the day
3. Follow the instructions above to set the day, hour and minutes.
4. Once you have set the minutes and pressed the TIME-button, the date and time are set and the time is displayed.
5. If you want to change the date and time, press and hold the TIME-button for approx. 3 seconds into the setting mode. Press TIME-button again and then the year number starts to flash. Now you can enter the new values as described above.

3. Checklist for Taking a Reliable Measurement

-  Avoid activity, eating or smoking immediately before the measurement.
-  Sit down on a back-supported chair and relax for 5 minutes. Keep the feet flat on the floor and do not cross your legs.
-  Always measure on the same wrist (normally left). It is recommended that doctors perform double arm measurements on a patient's first visit in order to determine which wrist to measure in the future. The wrist with the higher blood pressure should be measured.
-  Remove any items of clothing and your watch, for example, so that your wrist is free.

-  Always ensure that the cuff is positioned correctly, as shown in the pictures illustrated on the cuff.
-  Wrap the cuff on your wrist about 0.5 - 1 cm above your hand. Fit the cuff comfortably but not too tight. The cuff will cover a wrist circumference of 13.5 - 21.5 cm (5.25 - 8.5 inches).

4. Taking a Blood Pressure Measurement

1. Press the ON/OFF button (1) to start the measurement.
2. The cuff will now pump up automatically. Support your arm in a relaxed position with the device at the same height as your heart. Relax, do not move and do not tense your arm muscles until the measurement result is displayed. Breathe normally and do not talk.
3. The cuff fit check (18) on the display indicates that the cuff is perfectly placed. If the icon (18)-A appears, the cuff is fitted suboptimally, but it is still ok to measure.
4. When the correct pressure is reached, the pumping stops and the pressure falls gradually. If the required pressure was not reached, the device will automatically pump some more air into the cuff.
5. During the measurement, the pulse indicator (15) flashes in the display.
6. The result, comprising the systolic (12) and the diastolic (13) blood pressure and the pulse (14) is displayed. Note also the explanations on further display symbols in this booklet.
7. When the device has finished measuring, remove the cuff.
8. Switch off the device. (The monitor does switch off automatically after approx. 1 min.).

-  You can stop the measurement at anytime by pressing the ON/OFF button (e.g. if you feel uneasy or an unpleasant pressure sensation).

How not to store a reading

As soon as the reading is displayed press and hold the ON/OFF-button (1) until «M» (16) is flashing. Confirm to delete the reading by pressing the TIME-button (4).

-  «CL» is displayed when the reading is deleted from the memory successfully.

How do I evaluate my blood pressure?

The triangle on the left-hand edge of the display (10) points at the range within which the measured blood pressure value lies. The value is either within the optimum (green,) elevated (yellow) or high (red) range. The classification corresponds to the following ranges defined by International guidelines (ESH, AHA, JSH). Data in mmHg.

Range		Systolic	Diastolic	Recommendation
1	blood pressure too high	>135	>85	Seek medical advice
2	blood pressure elevated	130-135	80-85	Self-check
3	blood pressure optimum	<130	<80	Self-check

The higher value is the one that determines the evaluation. Example: a blood pressure value of 140/80 mmHg or a value of 130/90 mmHg indicates "blood pressure too high".

5. Appearance of the irregular heartbeat (IHB) symbol

This symbol ⑨ indicates that an irregular heartbeat was detected. In this case, the measured blood pressure may deviate from your actual blood pressure values. It is recommended to repeat the measurement.

Information for the doctor in case of repeated appearance of the IHB symbol:

This device is an oscillometric blood pressure monitor that also measures the pulse during blood pressure measurement and indicates when the heart rate is irregular.

6. Data Memory

This device automatically stores the last 30 measurement values.

Viewing the stored values

Press the M-button ③ briefly, when the device is switched off. The display first shows «M» ⑩ and «A», which stands for the average of all stored values. Pressing the M-button again displays the previous value. Pressing the M-button repeatedly enables you to move from one stored value to another.

-  Blood pressure readings with suboptimal cuff fit ⑮-A are not considered in the average value.
-  Pay attention that the maximum memory capacity of 30 memories is not exceeded. When the 30 memory is full, the oldest value is automatically overwritten with the 31th value. Values should be evaluated by a doctor before the memory capacity is reached - otherwise data will be lost.

Clearing all values

If you are sure that you want to permanently remove all stored values, hold down the M-button (the device must have been switched-off beforehand) until «CL ALL» appears and then release the button. To permanently clear the memory, press the TIME button while «CL ALL» is flashing. Individual values cannot be cleared.

Cancel deletion: press ON/OFF-button ① while «CL ALL» is flashing.

7. Battery Indicator and Battery change

Low battery

When the batteries are approximately $\frac{3}{4}$ empty the battery symbol ⑰ will flash as soon as the device is switched on (partly filled battery displayed). Although the device will continue to measure reliably, you should obtain replacement batteries.

Flat battery - replacement

When the batteries are flat the battery symbol ⑰ will flash as soon as the device is switched on (flat battery displayed). You cannot take any further measurements and must replace the batteries.

1. Open the battery compartment ⑥ at the back of the device.
 2. Replace the batteries - ensure correct polarity as shown by the symbols in the compartment.
 3. To set date and time, follow the procedure described in «Section 2.»
-  The memory is cleared and date and time must be reset - the year number therefore flashes automatically after the batteries are replaced.

Which batteries and which procedure?

-  Use two new, long-life 1.5V, size AAA alkaline batteries.
-  Do not use batteries beyond their date of expiry.
-  Remove batteries if the device is not going to be used for a prolonged period.

Using rechargeable batteries

You can also operate this device using rechargeable batteries.

-  Only use «NiMH» type reusable batteries.
-  Batteries must be removed and recharged when the flat battery symbol appears. They should not remain inside the device as they may become damaged (total discharge as a result of low use of the device, even when switched off).
-  Always remove the rechargeable batteries if you do not intend to use the device for a week or more.
-  Batteries cannot be charged in the blood pressure monitor. Recharge batteries in an external charger and observe the information regarding charging, care and durability.

8. Error Messages

If an error occurs during the measurement, the measurement is interrupted and an error message, e.g. «Err 3», is displayed.

Error	Description	Potential cause and remedy
«Err 1 » ⑧	Signal too weak	The pulse signals on the cuff are too weak. Re-position the cuff and repeat the measurement.
«Err 2» ⑱- B	Error signal	During the measurement, error signals were detected by the cuff, caused for instance by movement or muscle tension. Repeat the measurement, keeping your arm still.
«Err 3» ⑱- C	Abnormal cuff pressure	An adequate pressure cannot be generated in the cuff. A leak may have occurred. Check that the cuff is correctly connected and is not too loose. Replace the batteries if necessary. Repeat the measurement.
«Err 5»	Abnormal result	The measuring signals are inaccurate and no result can therefore be displayed. Read through the checklist for performing reliable measurements and then repeat the measurement.

«HI»	Pulse of cuff pressure too high	The pressure in the cuff is too high (over 299 mmHg) or the pulse is too high (over 200 beats per minute). Relax for 5 minutes and repeat the measurement.
«LO»	Pulse too low	The pulse is too low (less than 40 beats per minute). Repeat the measurement.

** Please immediately consult your doctor, if this or any other problem occurs repeatedly.*

9. Safety, Care, Accuracy Test and Disposal

Safety and protection

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device may only be used for the purposes described in these instructions. The manufacturer cannot be held liable for damage caused by incorrect application.
- This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in the «Technical Specifications» section.
- Protect it from:

-water and moisture

- extreme temperatures
- impact and dropping
- contamination and dust
- direct sunlight
- heat and cold

- The cuffs are sensitive and must be handled with care.
- Do not exchange or use any other kind of cuff or cuff connector for measuring with this device.
- Only pump up the cuff once fitted.
- Do not use this device close to strong electromagnetic fields such as mobile telephones or radio installations. Keep a minimum distance of 3.3 m from such devices when using this device.
- Do not use this device if you think it is damaged or notice anything unusual.
- Never open this device.
- If the device is not going to be used for a prolonged period the batteries should be removed.
- Read the additional safety information provided within the individual sections of this instruction manual.
- The measurement result given by this device is not a diagnosis. It is not replacing the need for the consultation of a physician, especially if not matching the patient's symptoms. Do not rely on the measurement result only, always consider other potentially occurring symptoms and the patient's feedback. Calling a doctor or an ambulance is advised if needed.

Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed. Be aware of the risk of strangulation in case this device is supplied with cables or tubes.

Device care

Clean the device only with a soft, dry cloth.

Cleaning the cuff

Carefully remove spots on the cuff with a damp cloth and soapsuds.

Accuracy test

We recommend this device is tested for accuracy every 2 years or after mechanical impact (e.g. being dropped). Please contact your local distributor to arrange the test (see foreword).

Disposal

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

10. Conditions for service under guarantee

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 3 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: electrode cups, elastic bands, power supply cables etc. nor regular checks and maintenance (calibration).
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@sanicare.it.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;
- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect

- use, installation or transport of the product;
- components or accessories subject to rapid wear and tear (electrode cups, elastic bands);
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and profit-making purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled

in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

e-mail: info@sanicare.it

www.sanicare.it

Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

11. Technical Specifications

Operating conditions:

10 -40 °C / 50-104 °F

15 - 90% relative maximum humidity

Storage conditions:

-20-+55°C/-4-+131°F

15 - 90% relative maximum humidity

Weight:

130g (including batteries)

Method:

oscillometric, corresponding to Korotkoff method: Phase I systolic, Phase V diastolic

Measurement range:

40 - 255 mmHg- blood pressure

40 - 199 beats per minute - pulse

Cuff pressure display range:

0- 299 mmHg

Resolution:

1 mmHg

Static accuracy:

pressure within ± 3 mmHg

Pulse accuracy:

$\pm 5\%$ of the readout value

Voltage source:

2x 1.5 V alkaline batteries; size AAA

Battery lifetime:

approx. 700 measurements (using new batteries)

IP Class:

IP22

Reference to standards:

EN1060-1 /-3/-4; IEC60601-1;

IEC 60601-1-2 (EMC); IEC 60601-1-11

Expected service life:

Device: 3 years or 10000 measurements

Accessories: 2 years

Accessories Included:

Instruction manual, 2 x AAA Batteries

This device complies with the requirements of the Medical Device Directive 93/42/EEC. Technical alterations reserved.

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SECURITY OF ALKALINE BATTERIES INTO PRODUCTS

1. Fire - Never delete the batteries throwing them in a fire. they could break and release the substances contained inside them.
2. If ingested or inserted, seek immediate medical assistance.
3. Do not store or leave the batteries in places accessible by children.
4. Remove old or weak batteries from appliances.
5. Do not mix old and new batteries in the same apparatus. There may be losses. The replacement of all batteries in a device must be done simultaneously.
6. Never mix different types of batteries in a single device. There may be losses.
7. Install carefully the batteries in the correct polarity (+ and -) according to the instructions provided by the manufacturer of the product
8. Do not leave the Products turn on with batteries when they are depleted. There may be losses.
9. Remove the batteries from any device that must remain unused for long periods of time
10. Avoid contact of the batteries with metal objects or other batteries because they may experience losses or they could break.
11. Do not recharge batteries because there may be some losses or they could break.
12. Not groped to disassemble the batteries
13. Store batteries in a cool and dry. Avoid sudden changes in temperature. Store the batteries in their original packaging before use.

14. RESPECT THE ENVIRONMENT - DISPOSE OF EXHAUSTED BATTERIES IN THE CONTAINERS

The battery and the possible Electronic Equipment -AEE- or anything that contains it must not be disposed as urban waste, but in compliance with Directive 2002/96/EC, 2006/66/EC and of National Laws transposition for the disposal of the product end of life, and must be subjected to a separate collection to avoid damage to the environment. Failure to comply with the above is sanctioned by Law.

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Via/Piazza/ Street/Square: _____ N°: _____

Città/ City and State: _____ CAP/ Postal Code: _____

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Microlife Corporation